eufy

For FAQs and more information, please visit:

[www.eufylife.com](http://www.eufylife.com)

Health & Wellness

## Owner's Manual

Smart Scale P2

T9148 51005003207 V01

## Table of Contents

02	Important Safety Instructions
07	What's In the Box
07	Inserting the Batteries
08	Selecting the Measuring Unit
11	Starting a Measurement
12	Weight-Only Measurement
13	Using the Measuring Tape
13	Maintenance
14	Troubleshooting
17	Specifications
18	Customer Service

English	01
Español	21
Deutsch	41
Français	61
Türkçe	79
Русский	97

## Important Safety Instructions

- This device is powered by 4 AAA batteries. Do not attempt to recharge non-rechargeable batteries as this can cause the batteries to leak or explode.
- Do not mix and use different types of batteries together, e.g. lithium-ion (Li-ion) and nickel-metal hydride (NiMH) batteries.
- Do not mix new and used batteries.
- If the device will not be used for an extended period of time, remove the batteries before storing the device in a cool and dry place.
- This device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of this appliance in a safe way and understand the hazards involved. Children should not play with the appliance. Cleaning and user maintenance should not be done by children without supervision.

02 EN

- Used batteries should be placed in a sealed plastic bag and disposed of safely according to local environmental regulations.

- Read all instructions before using the device. Keep the manual for future reference.
- Only use the device as described in this user manual. Any misuse may result in electric shocks, burns, start a fire, or cause other unexpected hazards.
- Be aware that this device is designed for individuals to measure their weight and to monitor their body composition. Any information provided by this device is in no way meant to treat, cure, or prevent diseases or illnesses from happening.
- DO NOT use this device if you are pregnant or if you suspect you may be pregnant. The effects of this device on a fetus are unknown and untested.
- DO NOT use this device if you use a wearable or implanted electronic device such as a defibrillator or pacemaker.
- This device should not be used by individuals who are acutely or chronically ill, suffering from diseases or taking medications that can affect water levels. The accuracy of the readings for these individuals has not been verified. Specific medical advice should be obtained from a physician.
- The device is equipped with a data transmission function that emits electromagnetic energy. Portable Radio Frequency (RF) communications equipment (including peripherals such as antenna cables and external antennas) should not be used within 30 cm (12 inches) of the device, otherwise, performance may be affected.
- Do not modify this device.
- Do not expose this device to extreme temperatures, direct sunlight, moisture, or corrosive environments.
- To prevent slipping, do not step on the device after bathing or showering.

EN 03

- This device should not be used when it is adjacent to or stacked on top of another device because it could cause the device to operate improperly. If such use is necessary, this device and the other device should be observed to verify that they are operating normally.

- If the batteries start to leak, put on your gloves and remove the leaking batteries from the battery door immediately. Place them in a clear plastic bag and recycle. If you don't wear gloves, you should wash your hands thoroughly with soap afterwards. In case of eye contact, seek medical attention immediately.

### Usage Tips:

To ensure measurement accuracy, follow the instructions below when you start your measurement:

- Place the scale on a flat, hard surface. Soft surfaces such as carpeting will affect the performance of the scale.
- Step barefoot onto the center of the scale. Stand still and keep full contact with the electrodes until the measurement is complete.
- For the most consistent measurements, weigh yourself at the same time of day and location.
- Avoid weighing yourself immediately after strenuous exercise, sitting in the sauna, taking a bath, or drinking and eating.
- The condition of the skin on the bottom of your feet can affect the measurement. The natural effects of aging or activity can make this skin hard. Take the reading with clean, slightly damp feet for best accuracy.
- Body fat percentage estimates will vary with the amount of water in the body, and can be affected by dehydration or overhydration due to factors such as alcohol consumption, menstruation, illness, intense exercise etc.
- Considering the vast individual differences among children under 16, the measurement of body fat or other body composition measurements may be less accurate. Their weight-only measurement should remain accurate.

## Notice

This symbol on the device or its packaging indicates: Do not dispose of electrical appliances as unsorted municipal waste; use separate collection facilities. Contact your local authority for information regarding the collection systems available. If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being. When replacing old appliances with new ones, the retailer is legally obliged to take back your old appliance for disposal at least free of charge.

## FCC Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**Warning:** Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

**Note:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: (1) Reorient or relocate the receiving antenna. (2) Increase the separation between the equipment and receiver. (3) Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. (4) Consult the dealer or an experienced radio / TV technician for help.

The following importer is the responsible party (for FCC matters only).

Company Name: POWER MOBILE LIFE, LLC

Address: 400 108th Ave NE Ste 400, Bellevue, WA 98004-5541

Telephone: +1 (800) 988 7973

### RF Exposure Compliance Statement

This equipment complies with FCC/IC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20cm between the radiator and your body.

### ISED Compliance Statement

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

### Declaration of Conformity

This product complies with the radio interference requirements of the European Community. Hereby, Anker Innovations Limited declares that the radio equipment type T9149 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <https://www.eufylife.com/>.

### GB Declaration of Conformity

Hereby, Anker Innovations Limited declares that the product type T9149 is in compliance with Radio Equipment Regulations 2017 and The Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012. The full text of the GB declaration of conformity is available at the following internet address:

<https://www.eufylife.com/>

The following importer is the responsible party (contract for EU and UK matters)

Anker Technology (UK) Ltd

Suite B, Fairgate House, 205 Kings Road, Tyseley, Birmingham, B11 2AA, United Kingdom

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
Georg-Mueller-Strasse 3, 80807 Munich, Germany

Anker Innovations Limited

Room 1318-19, Hollywood Plaza, 610 Nathan Road, Mongkok, Kowloon, Hong Kong

06 EN
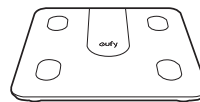
## What's In the Box

Smart Scale

AAA Battery (x4)

Soft Measuring Tape  
(60in / 150cm)


Owner's Manual &  
Other Documents

## Inserting the Batteries

- 1 Turn the scale over and place it on a flat surface.
- 2 Open the battery door and insert 4 AAA batteries. Make sure the positive and negative ends are facing the correct polarity directions marked in the battery compartment.  
- The digits will show on the LED display.
- 3 Press down to close the battery door. Make sure the battery door is fixed in place.

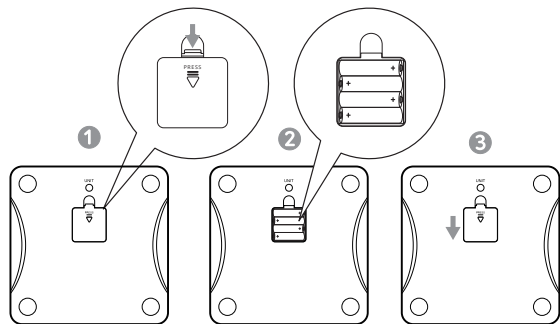EN 07


-  will appear on the LED display when the batteries are low. Replace all batteries at the same time. Do not mix new and used batteries.
- The scale turns off automatically after being idle for 10 seconds.

## Selecting the Measuring Unit

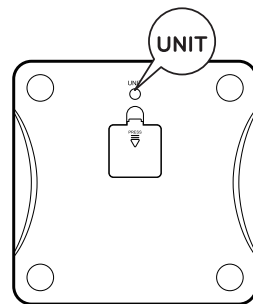Press the **UNIT** button at the bottom of the scale to cycle through and select the measuring unit: **lb** (Pound) / **kg** (Kilogram) / **st** (stone). The measuring unit can be selected and saved within the EufyLife app.

08 EN


- To restore factory settings, press and hold the **UNIT** button at the bottom of the scale for 10 seconds.

## Installing the EufyLife App and Pairing

- 1 Download the EufyLife app from the App Store (iOS devices) or Google Play (Android devices).
- 2 Enable Bluetooth on your smartphone or tablet.
- 3 Open the EufyLife app and follow the on-screen instructions to pair the scale with your device.

EN 09

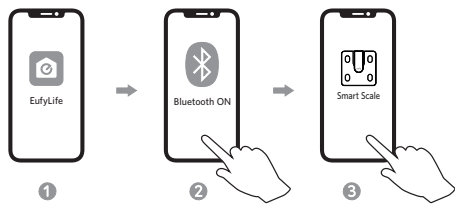- Once paired, the app and scale will automatically connect the next time the app is opened.

 - This scale is compatible with devices running iOS 10.0 and above or Android 5.0 and above.

- To complete the setup, you need to enable Bluetooth on your smartphone or tablet and allow the EufyLife app to use your location.

- Do not pair the scale from your smartphone's Bluetooth menu. Use the EufyLife app to pair the scale with your smartphone.


 Bluetooth Status Icon	Status Icon	Indicates
	Flashing	Pairing
	Solid	Connected

 Wi-Fi Status Icon	Status Icon	Indicates
	Off	Not connected
	Flashing slowly	Ready for connecting
	Flashing quickly	Connecting
	Solid	Connected

## Starting a Measurement

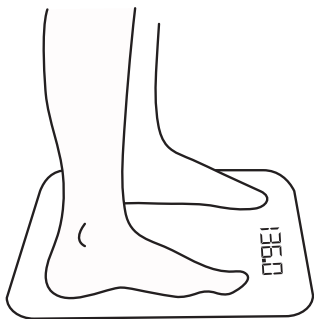Before taking your first measurement, you will need to create and set up a user ID in the EufyLife app.

- 1** Place the scale on a hard and flat surface.
- 2** Open the EufyLife app and select "Add Device" on the homepage. Follow the in-app instructions to add the smart scale.
- 3** Step barefoot onto the center of the scale. Stand still and keep full contact with the electrodes until the measurement is complete.
  - An analysis of your weight and your body's composition will begin automatically.
  - After measuring, your weight and other available information will be shown on the LED display. You can find detailed body composition analysis within the EufyLife app.


## Weight-Only Measurement

This scale can also be used as a conventional weight-reading scale.

1. Open the EufyLife app and select "Simple mode".
2. Step on the scale to measure your current weight.

## Using the Measuring Tape

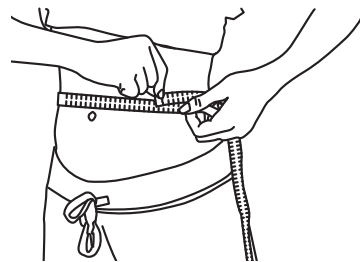With the measuring tape, you can get accurate measurements according to your requirements.

## Maintenance

Regularly clean the scale with a dry or slightly dampened soft cloth.

- This scale is waterproof, but it is not recommended to wash it with water for a long time or immerse it in water.
- Never use detergents to clean the scale as this may cause it to malfunction or discolor.

## Troubleshooting

### Error Message

Error Display	Description	Solution
	Overload. The device will power off.	Stop using this scale for measurement.
	Low battery. The device will power off.	Replace all batteries at the same time.

Error Display	Description	Solution
	Wi-Fi setup failure.	<ul style="list-style-type: none"><li>• Open the EufyLife app and follow the on-screen instructions to reconnect the scale with your device.</li><li>• Move the scale to a place near your router before the Wi-Fi setup.</li><li>• Press and hold the <b>UNIT</b> button for 10 seconds to restore factory settings, then restart the Wi-Fi connection.</li></ul>

### When Measuring

Issue	Cause	Solution
Abnormal measurement results: <ul style="list-style-type: none"><li>• Too high</li><li>• Too low</li><li>• Huge difference between two recent measurements</li></ul>	<ul style="list-style-type: none"><li>• Incorrect posture.</li></ul>	Step onto the scale barefoot and stand still.
	<ul style="list-style-type: none"><li>• Device is placed on soft carpeting or uneven surface.</li></ul>	Place the device on a hard and flat surface.

Issue	Cause	Solution
LED display does not turn on after stepping onto the scale.	• Batteries not installed.	• Install 4 new AAA batteries.
	• Worn batteries.	• Replace all batteries at the same time.
Scale does not proceed to analyse body composition.	• Taking a measurement while wearing shoes or socks.	• Step barefoot onto the center of the scale. Stand still and keep full contact with the electrodes until the measurement is complete.
	• Using the weight-only measuring mode in the EufyLife app.	• Change the measuring mode in the EufyLife app.

## Specifications

Model	T9148
Power	6V (4 × AAA batteries)
LED Display	2.5 × 1.4 in / 64 × 35 mm
Measurement Units	Pound / Kilogram / Stone
Unit Increments	0.1lbs / 0.05kg
Body Fat Unit Increments	0.1 %
Body Water Unit Increments	0.1 %
Muscle Mass Unit Increments	0.1 %
Measurement Range	0.2lbs - 397lbs / 0.1kg - 180kg / 0st: 0.2lbs - 28st: 5lbs
Activation Weight	3 kg
Working Environment	Temperature: 10 - 40 °C Humidity: <= 90% RH
Storage Environment	Temperature: -10 - 50 °C Humidity: <= 90% RH
Operating Frequency Range	2400 MHz - 2500 MHz
Waterproof Degree	IPX5
Transmitted Power	Wi-Fi: 18 dBm Bluetooth: 4 dBm

Bluetooth Version	5.0
Wi-Fi Version	2.4 G

## Customer Service

For tutorial videos, FAQs, manuals, and more information, please visit: [support.eufylife.com](https://support.eufylife.com).

### Contact Us

✉ [support@eufylife.com](mailto:support@eufylife.com)

💬 [www.eufylife.com](https://www.eufylife.com)

☎ (US) +1 (800) 988 7973 Mon-Fri 6:15 - 17:00 (PT);  
Sat-Sun 6:30-16:30 (PT)

(UK) +44 (0) 1604 936200 Mon-Fri 6:00 - 17:00;  
Sat-Sun 6:00-14:00

(DE) +49 (0) 69 9579 7960 Mo-Fr 8:00 - 16:00

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For more support phone numbers, please visit:

[support.eufylife.com/s/phonecontactus](https://support.eufylife.com/s/phonecontactus)

Or scan the QR code below:

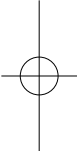
---

 @EufyOfficial  @EufyOfficial  @EufyOfficial  @EufyOfficial

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